

FIT - Specimen Collection Instructions

Do not collect sample during menstrual period or when suffering from hemorrhoids.
Pass urine **before** placing the collection paper on the water.

1




Unfold the collection paper and place on the water in the toilet bowl just before sitting on the toilet to collect the sample. The paper will float on the water for about 5 minutes.

2



Deposit stool sample on top of collection paper.

3

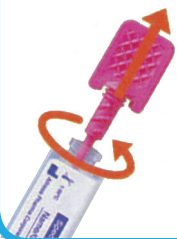


Unscrew the square tab and pull out the sample probe.

4



Pull sample probe across stool surface until the grooved portion of the probe has collected a sample of stool.

Sufficient sample is enough to cover the grooved portion of the sample probe.

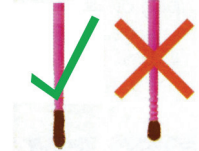

5



Close sample container(s) by inserting the sample probe and screwing the square tab tightly.

6



Fill in the collection date on the container label.

Place container(s) into the plastic bag provided. Return both containers to the nearest BC Biomedical Patient Service Centre within 7 days of collection.

7

After collection of sample, flush the paper with the stool

Important Precautions

- The sample containers are for single use only.
- The samples should be collected on 2 separate days or bowel movements.
- Store the sample containers in a cool dark place, before and after using
- Do not discard or spill liquid in sample container.
- Do not use the sample container for any other purpose than intended.
- This liquid contains a preservative. Keep the sample container out of reach of children and or pets.

FIT - 10.18.2011